
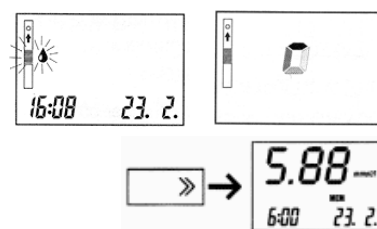


Glucotrend[®]2

Basic Meter Instructions

Performing a test

- Ensure you have set the time and date on your Meter before testing. Please refer to "Entering Settings" in the User's Guide.
- Turn the meter on by pressing centre button
- To check that all display elements are correct, wait for all three digits 888 to completely appear on the screen
- Check that the code on the screen matches the code number printed on the test strip vial. Remember to change the code calibration chip every time you open a new vial of test strips.
- Remove test strip from the vial. Replace vial cap.
- Insert the strip in the test slot in the direction shown by arrows until it clicks into place.
- The 3 digit code number will disappear.
- After approximately one second, a small blood drop symbol flashes on display indicating readiness to apply blood to test strip.
- After just a few seconds, you will hear a beep tone and a rotating box will indicate that the reading is in progress. Result will show in 15 seconds.
- Results will be saved automatically, if time and date are set correctly.

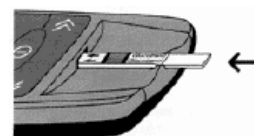


Coding the Glucotrend[®]2

- Remove the calibration coding chip from the test strip pack
- The coding chip slot is level with the bottom of the display window.
- Turn the chip so that the flat side with the code number inscribed on it, is facing upwards.
- Carefully insert the coding chip in the slot and press it in as far as it will go.
- If the monitor is off, switch it on.
- The bottom line of the display will show the time and date and a three digit code will be shown in the centre of the display. To the left of the screen, a test strip symbol and CODE will be flashing. A red dot flashes in the test-strip guide.
- Make sure that code on the monitor is the same as the code number on the pack of test strips.





Cleaning the Glucotrend® 2

- If your GLUCOTREND 2 is on, switch it off.
- Gently wipe across test strip guide using a damp non-linting cotton cloth to remove any traces of blood. Then using your thumb to lift it gently along its bottom edge, remove the test strip guide carefully.
- Clean the removed test strip guide with cold water and allow it to dry completely.
- Gently clean the exposed measuring window with a non-linting cotton cloth or cotton bud.
- Check carefully that test strip guide is completely dry, then re-insert.
- After cleaning, check the function of your meter using GLUCOTREND Plus Control G.

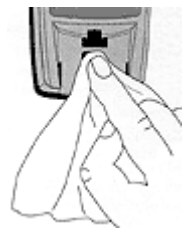

Changing the batteries

- Switch your GLUCOTREND 2 meter OFF.
- Gently press the battery compartment cover upwards with a thumb to open it, and then fold the cover on its hinges towards you. Do not try to remove battery cover completely.
- You can now see the twin battery holders. Insert one battery into each holder, with + upwards.
- Replace the battery compartment lid back over the batteries and press gently inwards to close it again.
- Press the ON key to check that the batteries are inserted properly.

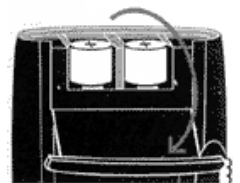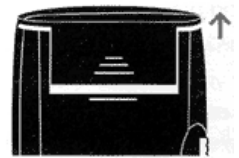

"For people with diabetes. Use only as directed and see your healthcare professional for medical advice."

Accu-Chek Integra is a trademark of Roche. © 2008 Roche Diagnostics



Roche Diagnostics Australia Pty. Ltd
ABN 29 003 001 205
Accu-Chek Enquiry Line 1800 251 816
www.accu-chek.com.au

ACCU-CHEK®
Live life. We'll fit in.